


# Injection Molding Small Complex Bioabsorbable Implants

Injection molding intricate components for medical device implantables is a difficult task on its own. Add absorbable polymer to the equation and it becomes many times more challenging. This article reviews several of the barriers to success in this area of component fabrication for implantable medical devices.

*By Tanner Hargens*

**B**ioabsorbable implants often pose several advantages to their non-absorbable counterparts. Absorbable materials do not require a second surgical procedure to remove the device, often reducing total charges and the rate of infection. These polymers may be tailored to serve an initial purpose after being implanted and then transfer responsibilities to the natural tissue in sync with the healing process. Drugs can also be impregnated within the polymer matrix to release at a rate controlled by the degradation of the material.

Bioabsorbables are widely utilized for cardiovascular applications. Examples include ligation and closure of blood vessels along with

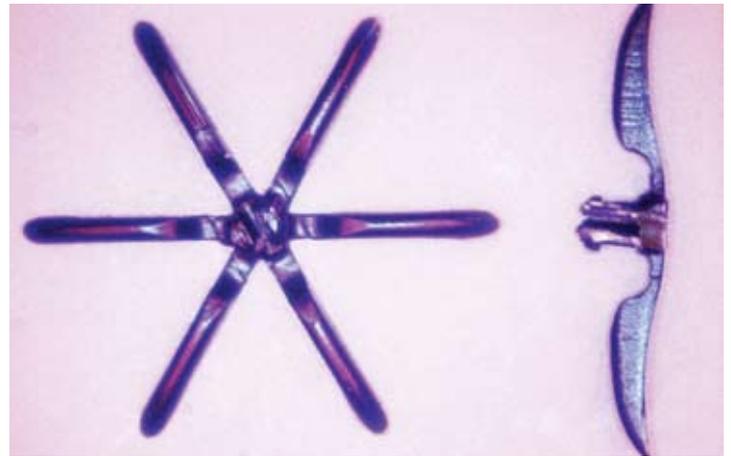

**Bioabsorbable implant with snap feature and six living hinges (part diameter: 14 mm)**

absorbable drug releasing stents and grafts. Many of these devices have small, complex, 3D structures and dynamic mechanical requirements (e.g., living hinges) that can only be formed by injection molding.

Injection molding first transitions a polymer from a solid state to a semi-liquid state by heating the material above a threshold temperature. With a reduced viscosity upon turning into a liquid, the polymer can be injected (or forced) into the cavity of a mold. The process must then be reversed (i.e., the material must cool and return to a solid state) before it can be removed from the mold.

Bioabsorbable polymers can be injection molded with extremely small features and thin sections; however, this requires special equipment and consideration. Absorbable polymer should be processed at the lowest temperatures possible with the least amount of time exposed to the elevated temperatures. Each bioabsorbable will have a different temperature threshold that, when surpassed, will reverse the polymerization process resulting in monomer formation. Introducing an extra monomer into the material matrix can immediately change the mechanical properties, alter the degradation kinetics, and even make the material toxic. Conventional injection molding equipment and processes are unable to adequately perform such controlled and rapid thermal cycling.

In addition, presses need much higher injection speeds and pressures to push melted polymer through tiny nozzles and flow channels. Higher injection speeds ensure the material will fill the part before it cools. Accurate screw and barrel design and injection control are essential. Molding very small parts generally is not recommended on machines



## Projected Capacitive Touch Display Provider




**If you need further assistance, please contact us.**

**DATA IMAGE CORPORATION**

www.dataimagelcd.com  
 heldred@dataimagelcd.com  
 office: 407-478-4064

Sizes Including		
Size		Resolution
3.5"	QVGA	320x240
4.3"	(MVA)WVGA	480x800
5"	WVGA	800x480
5.7"	VGA	640x480
7"	WVGA	800x480
7"	(IPS)WSVGA	1024x600
10.1"	WSVGA	1024x600
10.1"	WXGA	1280x800



### Injection molded bioabsorbable polymer implant

larger than 20 tons because it is difficult to control such small shot sizes, and lab machines, though small, are not usually designed for precision control of volume production.

Those that are capable of molding such small parts commonly utilize sprue and runner systems consuming 75% or more of the volume needed to produce one part. Medical grades of bioabsorbable resins can sell for upwards of \$13 per gram, so the losses must be tightly controlled to minimize the effect on the bottom line.

Prior to and after melt processing, these absorbable polymers should be kept dry and shielded from any humidity to extend material shelf life. This can be done by utilizing foil pouches that have had oxygen removed, desiccants added, and/or freezing temperatures administered to minimize the hydrolytic degradation. Every aspect of the manufacturing process must be closely monitored and proper procedures strictly followed. Due to the special attention that is needed with bioabsorbable parts, it is best to use an injection molding unit dedicated for absorbables with dedicated processes (e.g., cleaning with abrasive plastic cloth, dryer/freezer for storage, vacuum storage of molded parts, personnel understanding critical effects, etc.).

For these reasons, specialized systems have been developed for injection molding small complex bioabsorbable parts. Medical Murray ([www.medicalmurray.com](http://www.medicalmurray.com)) uses a molding machine called the Sesame, which was designed and patented (U.S. No. 6,267,580) by the company for creating small complex medical devices. The machine is designed to provide fast controlled injection and short material residence time. It uses a pneumatically driven vertical plunger for plasticating and injects with a tiny horizontal plunger (1.0- to 2.5-mm diameter) driven by a linear servo motor. Total injection time can be as short as 0.02 seconds. The Sesame's features include very low melted volume and controlled high speed

and high pressure injection. This allows the molding of sub-micro size parts with complex geometric features from materials such as bioabsorbable polymers that would otherwise degrade in standard equipment. Parts have been produced on the Sesame measuring just 0.1 mm and weighing only 2.5 mg, or one-twentieth of a single pellet.

*Tanner Hargens is a senior biomedical engineer and project leader for Medical Murray Inc. He is responsible for medical device engineering research and development, and specializes in implants, biomaterials, and minimally invasive procedures. Hargens also provides expertise related to commercializing medical products.*

### Ask the Experts!

Question about this article? Send it to [sean.fenske@advantagemedia.com](mailto:sean.fenske@advantagemedia.com) and the response will be shared in MDT's blog online.

## The World's Leading Manufacturer of Implantable Crystal Devices

Over 20 Years Experience in Medical Implantable and Non-Implantable Applications

**Medical Applications:** Pacemaker, Defibrillator, CRM, Neurostimulation, Telemetry, Insulin Pump, Glucose Meter

**Production:** Dedicated Medical Production Lines and Experienced in High Volume Consumer Products

**SWISS MADE QUALITY**

OV-7604-C7	CM8V	CC7A	CC7V/CM7V	CC5V	CC4V
32.768 kHz	32.768 kHz	24 to 30 MHz	32 to 200 kHz	32 to 200 kHz	30 to 1000 kHz
3.2x1.5x1.0	2.0x1.2x0.6	3.2x1.5x0.9	3.2x1.5x0.9	4.1x1.5x0.9	5.0x1.9x0.9
Oscillator			3.2x1.5x0.65		



MICRO CRYSTAL SWITZERLAND

Micro Crystal AG  
Muehlestrasse 14 CH-2540 Grenchen Switzerland  
Phones International: +41 32 655 82 82 US: +1 847 818 98 25  
[sales@microcrystal.com](mailto:sales@microcrystal.com) [www.microcrystal.com](http://www.microcrystal.com)

A COMPANY OF THE SWATCH GROUP